



Our Roomer S overhead lift is a portable lift that offers the possibility to transfer a patient from one room to another in one single motion without making modifications to doorway headers.

A Comfortable Experience

The Roomer S offers optimum comfort when used with our sling portfolio. The room-to-room function reduces the number of lifts required and minimizes strain on both the patient and the caregiver. It offers a soft-start, soft-stop and constant lifting speed regardless of the patient weight and direction, making lifting a comfortable experience. The lift has a weight capacity of 220kg/484lb and can be used with either a 2, 3 or 4-point suspension.

Unique Design

The Roomer S is suitable for rooms with sliding doors or fire doors as no modification to the doorway is required, thus making installation

much easier. For facilities, its unique design eliminates the need for turntables and rail switches. The Roomer S also features a lightweight quick release hanger bar in two different widths to accommodate most lifting situations. Charging can be done through the hand control for optimal efficiency.

Lift Intelligence

The Roomer S features an advanced technology platform that allows tracking of lift usage and visualization of collected data. This integrated lift intelligence allows data to be retrieved using our Data and Service Application on a tablet or computer. Detailed information can be viewed directly in the application or exported to other media for additional analysis.

HUMAN CARE SWEDEN (HQ)

Årstaängsvägen 21B
117 43 Stockholm
Phone: +46 8 665 35 00
Fax: +46 8 665 35 10
info.se@humancaregroup.com

HUMAN CARE AUSTRALIA

Unit 4, 49 Bertie Street
Port Melbourne, VIC 3207
Phone: +61 1300 132 736
Fax: +61 (03) 9645 4954
info.au@humancaregroup.com

HUMAN CARE CANADA

10-155 Colonnade Road
Ottawa, ON K2E 7K1
Phone: +613.723.6734
Fax: +613.723.1058
info.ca@humancaregroup.com

HUMAN CARE NEDERLAND

Elspeterweg 124
8076 PA Vierhouten
Phone: +31 577 412 171
Fax: +31 577 412 170
info.nl@humancaregroup.com

HUMAN CARE UNITED STATES

8006 Cameron Road, Suite K
Austin, TX 78754
Phone: +512.476.7199
Fax: +512.476.7190
info.us@humancaregroup.com

Accessories



Hanger bar

Quick release for optimal convenience. Order separately. Part no.: 55751 [45cm/17.7"] 55764 [55cm/21.7"]



Upholstered Cover

Protects user from hanger bar impact. Hanger bar sold separately. Part no. 50866



4-point suspension

Enables 4-point lifting. Placed on leg hook. Part no.: 55619



Cart

Simplifies transport of lift between rooms. Convenient sling storage. Part no.: 50346



Powered spreader bar

For use with clip slings in hospital or other care environments. Part no.: Small 92441, Medium 92442, Large 92443



Lift charger

Charger with holder. Easily attaches to wall. Part no. 50899-US



Distance Strap

Reinforced nylon thread and metal hooks. Extend lift strap up to 100cm/39.3". Part no. 20124-xx [specify length]




Scale

Medical grade scale. Max. capacity 320kg/705.5lb. ISO and EN compliant. Part no.: 50348, 50345 [Canada]



Handset

Two speed with up/down mode. Input for charging and USB interface. IPX5 water protection rating. Part no.: 71457

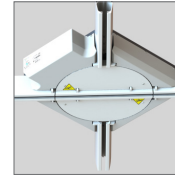


Gates and Turntables

Offering additional lift transfer options on rails. Part no. 55721 [Gate], 91884 [Turntable]



Secondary Control

Secondary control panel on back of lift for emergency stop function. Integrated feature.

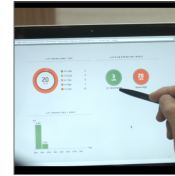

Insights

Lift intelligence that facilitates cloud based storage of collected data. Integrated feature.

Technical Specifications

Description	Roamer S
Part Number (Lift)	55300H-[Country Code]
Part Number (Hanger Bar)	55751 [45cm/17.7"] 55764 [55cm/21.7"]
Dimensions (W x H x L)	40cm x 16cm x 30cm 15.7" x 6.2" x 11.8"
Safe Working Load	220kg/484lb
Product Weight (Lift)	9.9kg / 21.7lb
Product Weight (Hanger Bar)	1.3kg/2.8lb [55751] or 1.5kg/3.3lb [55764]
Battery	2.3Ah, 2x12V
Battery Charger	Input: 100-240V, AC 50-60 Hz Output: 27.8V, 0.8A DC
IP Class	Lift: IPX4 Hand Control: IPX5
Lifting Distance and Speed	Distance: 208cm/81.9" Speed: 1.8m/min - 3m/min — 70"/min - 118"/min
Approvals	ROHS-2 approved, Tested and approved according to ISO 10535:2006 ETL Listed
Emergency Lowering	Electronic and Mechanical

This product is manufactured in accordance with all relevant ISO standards.



Learn more about this product



Country Codes: AU-Australia, GB-Great Britain, EU-Europe, JP-Japan, NA-North America

Follow Us on Social Media



facebook.com/humancareworld



instagram.com/humancaregroup



linkedin.com/company/human-care-ab



pinterest.com/humancaregroup



youtube.com/channel/UC7fjINh3zTjRxOP1IDf2Da0