

1700



The 1700 Convertible Chair is a no-lift, multi-functional device that helps caregivers to more efficiently and effectively perform safe patient handling practices that can help to improve patient satisfaction, safety, dignity, and reduce pain. Its Tilt-n-Space capability and therapeutic cushioning allows increased blood circulation helping to prevent skin breakdown and other complications of immobility.

Unlimited Repositioning

The 1700 Barton Convertible Chair is a portable, efficient and versatile transfer and positioning solution for bariatric patients that helps move a patient from a supine to a seated position easily and safely. By simply getting the patient out of the bed faster, it promotes faster recovery and can help reduce the length of stay in both higher acuity departments as well as the overall length of stay in a hospital. A non-lifting alternative, it is an ideal option for patients with pressure injuries, pain or those who are simply afraid of being lifted. The caregiver can continuously and frequently reposition a patient to increase tolerance for long-term seating with the recline or Tilt-n-Space functions.

Easy Transfers

The 1700 Barton Convertible chair when

used with the patient transfer system allows supine transfers in four easy steps (see illustrations). The ability to use the 1700 improves caregiver safety while eliminating unnecessary patient transfers in addition to allowing for transport to physical therapy and other ancillary services

Quality

Our products are designed and manufactured to meet all the requirements of relevant standards. At every stage of the development process, every single product is given extensive attention in design, development and manufacturing. This means that a user of this product feels confident with its reliability and the caregiver feels confident with its safety.

HUMAN CARE SWEDEN (HQ)

Årstaängsvägen 21B
117 43 Stockholm
Phone: +46 8 665 35 00
Fax: +46 8 665 35 10
info.se@humancaregroup.com

HUMAN CARE AUSTRALIA

Unit 4, 49 Bertie Street
Port Melbourne, VIC 3207
Phone: +61 1300 132 736
Fax: +61 (03) 9645 4954
info.au@humancaregroup.com

HUMAN CARE CANADA

10-155 Colonnade Road
Ottawa, ON K2E 7K1
Phone: +613.723.6734
Fax: +613.723.1058
info.ca@humancaregroup.com

HUMAN CARE NEDERLAND

Elspeterweg 124
8076 PA Vierhouten
Phone: +31 577 412 171
Fax: +31 577 412 170
info.nl@humancaregroup.com

HUMAN CARE UNITED STATES

8006 Cameron Road, Suite K
Austin, TX 78754
Phone: +512.476.7199
Fax: +512.476.7190
info.us@humancaregroup.com

Accessories



PTS System

Patient transfer system used for transfer to/from a bed to/from the chair. Part no. b2013a




Transfer Sheets

Standard or bariatric sizing. Part no. b1410 (standard), b1420 (bariatric).



Castors

Front: 3 Position, Rear: 2 Position. Replacement Part no.: l-b10x2-8049 (front), l-b10x2-8050 (rear)

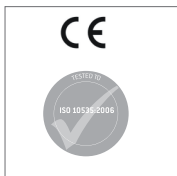


Manufactured to all relevant standards, including

Transfers




Step 1

Place the chair by the bed in the stretcher position.

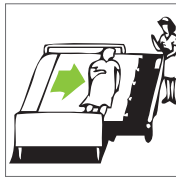


Step 2

Attach the bed sheet to the transfer bar.

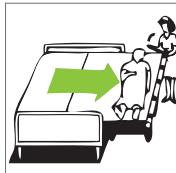

Step 3

Transfer the patient by turning the handle.

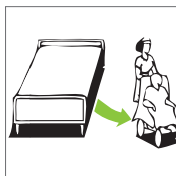

Step 4

Incline the chair to transport patient.

Features



Headrest

Adjustable in height, depth and angle for proper support.



Tilt-n-Space

Designed for pressure relief and positioning. Moves position without affecting body angle of patient.



Adjustable Leg/Footrest

Adjustable to 90° angle and 6" for added length and comfort.



Anti-Bacterial Fabric

Anti-bacterial, anti-microbial, breathable material that can be easily cleaned or disinfected in seconds.

Technical Specifications

1700			
Part No.	B1071-TS-GEL	Overall Height (Prone)	70cm / 27.5"
Seat Width	65cm / 25.5"	Overall Height (Wheelchair)	130cm - 142cm / 51" - 56"
Overall Width	89cm / 35"	Back Angle	0° to 70°
Overall Length (Prone)	184cm - 199cm / 72.5" - 78.5"	Tilting Angle	0° to 17°
Overall Length (Wheelchair)	115cm - 159cm / 45.5" - 62.5"	Castors	15cm / 6"
Headrest (Adjustable Height)	0cm - 12.7cm / 0" - 5"	Product Weight	93.4kg / 206lb
Footrest (Adjustable Height)	0cm - 12.7cm / 0" - 5"	Weight Capacity	318kg / 700lb



Watch our chair transfer video

This product is manufactured in accordance with all relevant ISO standards.



Learn More About This Product



Follow Us on Social Media



facebook.com/humancareworld



instagram.com/humancaregroup



linkedin.com/company/human-care-ab



pinterest.com/humancaregroup



youtube.com/channel/UC7fjInh3zTjRx0P1IDf2Da0

humancaregroup.com