
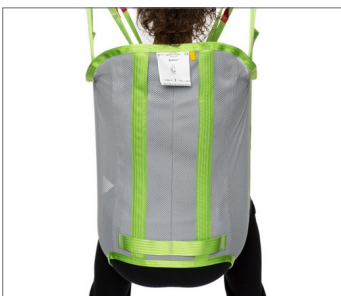


# bath sling

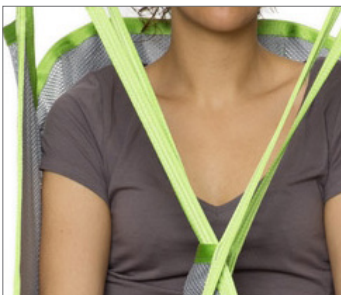

The Bath Sling can be used in most lifting situations and fits most patients but is most suitable for wet environments. The Polyester Net material quickly dries after use making it ideal for bathing, showering and rehabilitation in pool areas. The Bath Sling designed to give the patient extra support and is a good solution for those users with poor trunk stability.

#### Quick Dry Material

The Bath Sling is made of a Polyester Net fabric, a soft and thin material that breathes and dries quickly. It is ideal choice for bathing and showering, but can also be used in rehabilitation pools. The sling design creates an upright seating that ensures proper positioning when the patient is placed in a bath chair. Because the material has some elasticity, it can be left behind the patient for a short period of time.

#### Extra Support

The Bath Sling features a reinforced back and leg parts that offer additional support throughout the back area for those users with poor trunk stability. The Bath Sling can

be used with or without an optional waist belt that is available in an assortment of sizes with a Velcro or buckle closure.

#### Easy Care and Easy Fit

Our slings are manufactured to the highest quality standards. They feature reinforced straps and high strength nylon webbing around the sling to make them strong enough to resist even the most active use. Each sling label is colour coded to ensure quick and easy selection of the correct sling and the straps are colour coded for easy determination of required positioning.

#### HUMAN CARE SWEDEN (HQ)

Årstaängsvägen 21B  
117 43 Stockholm  
Phone:- +46 8 665 35 00  
Fax:- +46 8 665 35 10  
info.se@humancaregroup.com

#### HUMAN CARE AUSTRALIA

Unit 4, 49 Bertie Street  
Port Melbourne, VIC 3207  
Phone : +61 1300 132 736  
Fax : + 61 (03) 9645 4954  
info.au@humancaregroup.com

#### HUMAN CARE CANADA

10-155 Colonnade Road  
Ottawa, ON K2E 7K1  
Phone: +613.723.6734  
Fax: +613.723.1058  
info.ca@humancaregroup.com

#### HUMAN CARE NEDERLAND

Elspeterweg 124  
8076 PA Vierhouten  
Phone: +31 577 412 171  
Fax: +31 577 412 170  
info.nl@humancaregroup.com

#### HUMAN CARE UNITED STATES

8006 Cameron Road, Suite K  
Austin, TX 78754  
Phone: +512.476.7199  
Fax: +512.476.7190  
info.us@humancaregroup.com

## Accessories



### Extension Loops

Used to extend lift straps in order to alter lifting position or to hold leg parts together. Assorted lengths. Part no. 21180.\*



### Leg Padding and Supports

Ideal for users with sensitive skin or pain. Part no.: 21146-S, M, L Reinforced supports. One size. Part no. 25102



### Waist Belts

Velcro or buckle enclosure. Available in Junior and Adult sizing. Velcro Part no. 25145-\*\*\* Buckle Part no. 25146-\*\*\*



### Hanger Bars for 2-, 3- or 4-point Suspension

Part no. 57551, 55764 [2- or 3-Point], 55619 [4-Point - Attaches to leg hook]

## Lifting Situations



### Use to/from a bed

Suitable for transfers and repositioning to and from a bed. Can be left behind patient for short period of time.



### Use to/from a chair


Suitable for transfers and repositioning to and from a chair. Can be left behind patient for short period of time.



### Use for showering

Only Polyester Net fabric can be used for this application.

## Care and Features

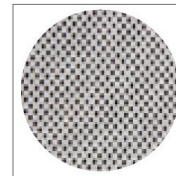


### Polyest Net

Soft, thin material that breathes and dries quickly. Ideal for bathing or showering.



### Care

Care instructions located on sling label. Use a hospital grade laundry bag while laundering.

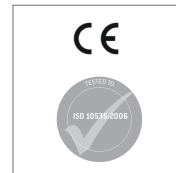


Tested to all relevant standards, including ISO 10535, hoists for the transfer of disabled persons.

This product is manufactured in accordance with all relevant ISO standards.



## Technical Specifications

Part Number	Size	Colour	A	B	C
25035-JS	Junior Small	JS	22cm / 8.7"	56cm / 22"	53cm / 22.9"
25035-JM	Junior Medium	JM	28cm / 11"	63cm / 24.8"	58cm / 22.8"
25035-S	Small	S	35cm / 13.8"	76cm / 29.9"	62cm / 24.4"
25035-M	Medium	M	40cm / 15.7"	86cm / 33.9"	62cm / 24.4"
25035-L	Large	L	44cm / 17.3"	98cm / 38.6"	63cm / 24.8"
25035-XL	Extra Large	XL	48 cm / 18.9"	102cm / 40.2"	64cm / 25.2"
Maximum User Weight	250kg/551lb				

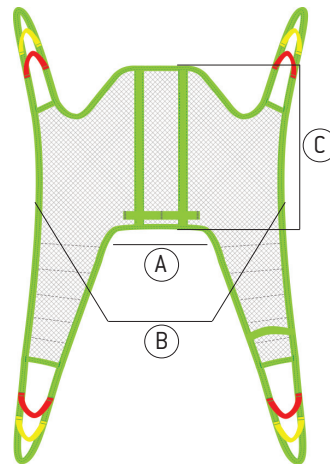

\* Available in assorted lengths: 10cm/3.9" [21180-10], 15cm/5.9" [21180-15], 20cm/7.8" [21180-20], 25cm/9.8" [21180-25], 35cm/13.7" [21180-35]

\*\* Junior Sizes: 40cm/15.7" [25145-40B], 50cm/19.7" [25145-50B], 60cm/23.6" [25145-60B]; Adult Sizes: 60cm/23.6" [25145-60], 70cm/27.6" [25145-70], 80cm/31.5" [25145-80], 90cm/35.4" [25145-90], 110cm/43.3" [25145-110], 120cm/47.3" [25145-120]

\*\*\* Junior Sizes: 40cm/15.7" [25146-40B], 50cm/19.7" [25146-50B], 60cm/23.6" [25146-60B]; Adult Sizes: 60cm/23.6" [25146-60], 70cm/27.6" [25146-70], 80cm/31.5" [25146-80], 90cm/35.4" [25146-90], 100cm/39.4" [25146-100], 110cm/43.3" [25146-110], 120cm/47.3" [25146-120]

Learn More About This Product



Follow Us on Social Media



facebook.com/humancareworld



instagram.com/humancaregroup



linkedin.com/company/human-care-ab



pinterest.com/humancaregroup



youtube.com/channel/UC7fjInh3zTjRxOP1IDf2Da0

humancaregroup.com