

# toileting slings



Our Toileting Slings feature a wide, divided leg partition that enables fast and simple application for dressing and undressing during toileting or bathing. They are suitable for patients with poor trunk stability and feature padding under the arms for extra comfort and are available with extra head support or without.

### Three Options

These Toileting slings are available in two low back models and one high back model with adjustable neck support for increased head stability. The wide partition of the leg parts makes the sling application fast and simple and enables dressing and undressing during toileting or bathing. All models can be used with or without an optional waist belt that is available in an assortment of sizes with a Velcro or buckle closure.

### Optimal Positioning

The low back models feature an upright seated position that makes transferring to and from a toilet much easier. The high back model features a reclined (backward tilted)

position that prevents the user from sliding out of the sling and is especially beneficial for those hypertonic users that have a tendency to overstress themselves (spastic performance) or those hypotonic users with lesser body control, dizziness or users with extended time spent in bed.

### Easy Care and Easy Fit

Our slings are manufactured to the highest quality standards. They feature reinforced straps and high strength nylon webbing around the sling to make them strong enough to resist even the most active use. Each sling label is colour coded to ensure quick and easy selection of the correct sling and the straps are colour coded for easy determination of required positioning.

#### HUMAN CARE SWEDEN (HQ)

Årstaängsvägen 21B  
117 43 Stockholm  
Phone: +46 8 665 35 00  
Fax: +46 8 665 35 10  
info.se@humancaregroup.com

#### HUMAN CARE AUSTRALIA

Unit 4, 49 Bertie Street  
Port Melbourne, VIC 3207  
Phone : +61 1300 132 736  
Fax : + 61 (03) 9645 4954  
info.au@humancaregroup.com

#### HUMAN CARE CANADA

10-155 Colonnade Road  
Ottawa, ON K2E 7K1  
Phone : +613.723.6734  
Fax: +613.723.1058  
info.ca@humancaregroup.com

#### HUMAN CARE NEDERLAND

Elspeterweg 124  
8076 PA Vierhouten  
Phone: +31 577 412 171  
Fax: +31 577 412 170  
info.nl@humancaregroup.com

#### HUMAN CARE UNITED STATES

8006 Cameron Road, Suite K  
Austin, TX 78754  
Phone: +512.476.7199  
Fax: +512.476.7190  
info.us@humancaregroup.com

## Accessories



### Extension Loops

Used to extend lift straps in order to alter lifting position or to hold leg parts together. Assorted lengths. Part no. 21180.\*



### Waist Belts

Velcro or buckle enclosure. Available in Junior and Adult sizing. Velcro Part no. 25145-\*\* Buckle Part no. 25146-\*\*\*



### Upholstered Cover

Protects user from hanger bar impact. Not for use in a wet environment. Part no. 50866



### Hanger Bars for 2-, 3- or 4-point Suspension

Part no. 57551, 55764 [2- or 3-Point - Quick Release], 55619 [4-Point - Attach to leg hook]

## Lifting Situations



### Use to/from a bed

Suitable for transfers and repositioning to and from a bed. Can be left behind patient for short period of time.




### Use to/from a chair

Suitable for transfers and repositioning to and from a chair. Can be left behind patient for short period of time.



### Use for showering


Only Polyester Net fabric can be used for this application.

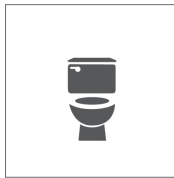

### Use for toileting


Only Polyester Net fabric can be used for this application.

## Care and Features

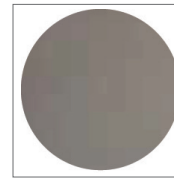


### Polyester

Durable, low-friction, easy to clean, apply and remove.

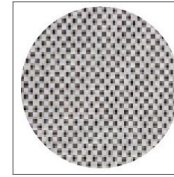


### Polyest Net

Soft, thin material that breathes and dries quickly. Ideal for bathing or showering.

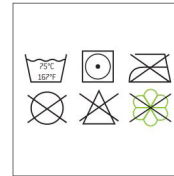

### Care

Care instructions located on sling label. Use a hospital grade laundry bag while laundering.

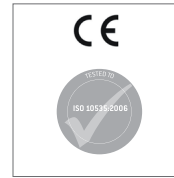

Tested to all relevant standards, including ISO 10535, hoists for the transfer of disabled persons.

## Technical Specifications

Part Number	Size	Colour	Model #25080/25085			Model #25090			Model #25080/25085	Model #25090
			A	B	C	A	B	C		
<b>Polyester Net</b>	<b>Polyester</b>									
-	25080-JS	Junior Small	JS	28cm / 11"	62cm / 24.4"	22cm / 8.7"	-	-	-	
25085-JM	25080-JM	Junior Medium	JM	28cm / 11"	67cm / 26.3"	26cm / 10.2"	28cm / 11"	73cm / 28.7"	70cm / 27.6"	
25085-S	25080-S	Small	S	32cm / 12.6"	74cm / 29.1"	30cm / 11.8"	32cm / 12.6"	82cm / 32.3"	74cm / 29.1"	
25085-M	25080-M	Medium	M	48cm / 18.9"	90cm / 35.4"	33cm / 13"	48cm / 18.9"	92cm / 36.2"	83cm / 32.7"	
25085-L	25080-L	Large	L	49cm / 19.3"	97cm / 38.2"	33cm / 13"	49cm / 19.3"	99cm / 39"	90cm / 35.4"	
25085-XL	25080-XL	Extra Large	XL	56cm / 22"	100cm / 39.4"	38cm / 15"	56cm / 22"	101cm / 39.8"	100cm / 39.4"	
Maximum User Weight	250kg/551lb									

\* Available in assorted lengths: 10cm/3.9" (21180-10), 15cm/5.9" (21180-15), 20cm/7.8" (21180-20), 25cm/9.8" (21180-25), 35cm/13.7" (21180-35)  
 \*\* Junior Sizes: 40cm/15.7" (25145-40B), 50cm/19.7" (25145-50B), 60cm/23.6" (25145-60B); Adult Sizes: 60cm/23.6" (25145-60), 70cm/27.6" (25145-70), 80cm/31.5" (25145-80), 90cm/35.4" (25145-90), 110cm/43.3" (25145-110), 120cm/47.3" (25145-120)  
 \*\*\* Junior Sizes: 40cm/15.7" (25146-40B), 50cm/19.7" (25146-50B), 60cm/23.6" (25146-60B); Adult Sizes: 60cm/23.6" (25146-60), 70cm/27.6" (25146-70), 80cm/31.5" (25146-80), 90cm/35.4" (25146-90), 100cm/39.4" (25146-100), 110cm/43.3" (25146-110), 120cm/47.3" (25146-120)

This product is manufactured in accordance with all relevant ISO standards.



Learn More About This Product



### Follow Us on Social Media



facebook.com/humancareworld



instagram.com/humancaregroup



linkedin.com/company/human-care-ab



pinterest.com/humancaregroup



youtube.com/channel/UC7fjInh3zTjR0P1Df2Da0