



Our lift accessories are manufactured to the highest quality standards and are designed to complement and increase the functionality of our lift systems. Simple to install and easy to use, these accessories help to adapt the lift solution to each individual need.

#### Hanger Bars

We offer multiple options in hanger bars including our newest 2-point quick release series - a lightweight yet durable hanger bar that connects in seconds! They feature an easy-to-read label that indicates the specifications, making it easy to choose the correct hanger bar. With extra large hooks and a high weight capacity, these hanger bars meet most lifting requirements. In addition, we offer 4-point systems to accommodate more robust patients and hanger bars specifically for fixed lift systems.

#### Medical Grade Scale

Our scale is approved for weighing patients in a medical environment. The scale is quick, accurate and easy to use and compliant with

all relevant standards.

#### 4-Point Powered Spreader Bar

This spreader bar is the optimal choice for hospitals and other health care environments who use slings that are affixed with clips. Controlled by the lift handset, the electric movement is smooth and comfortable for both the client and the caregiver.

#### Ceiling Lift Trolley

Designed to carry our portable lifts from one room or department to the next, our trolley features a non-slip surface, two directional and two rear locking casters. Easy to assemble it comes complete with a deep basket to carry slings and other lift accessories.

#### HUMAN CARE SWEDEN (HQ)

Årstaängsvägen 21B  
117 43 Stockholm  
Phone: +46 8 665 35 00  
Fax: +46 8 665 35 10  
info.se@humancaregroup.com

#### HUMAN CARE AUSTRALIA

Unit 4, 49 Bertie Street  
Port Melbourne, VIC 3207  
Phone : +61 1300 132 736  
Fax : + 61 (03) 9645 4954  
info.au@humancaregroup.com

#### HUMAN CARE CANADA

10-155 Colonnade Road  
Ottawa, ON K2E 7K1  
Phone : +613.723.6734  
Fax : +613.723.1058  
info.ca@humancaregroup.com

#### HUMAN CARE NEDERLAND

Elspeterweg 124  
8076 PA Vierhouten  
Phone: +31 577 412 171  
Fax: +31 577 412 170  
info.nl@humancaregroup.com

#### HUMAN CARE UNITED STATES

8006 Cameron Road, Suite K  
Austin, TX 78754  
Phone: +512.476.7199  
Fax: +512.476.7190  
info.us@humancaregroup.com

## Accessories



### Hanger Bars

Quick release series.  
Altair/Roomer S Part no. 55751, 55764. HeliQ/HeliQ PT Part no. 72754, 72755, 72757



### Hanger Bar - Fixed

Available in 4 sizes.  
Part no. 55611, 55612, 55613 and 54043



### Upholstered Cover

Protects user from hanger bar impact. Hanger bar sold separately. Part no. 50866



### 4-Point Suspension

Enables 4-point lifting on HeliQ/HeliQ PT. Part no. 70800



### 4-Point Suspension

Enables 4-point lifting on Altair/Roomer S. Placed on leg hook. Part no. 55619



### Scale

Medical grade scale can be used on all Human Care lifts. Part no. 50590 (Altair/Roomer S), 50591 (HeliQ), 50345 (NA only)



### Trolley

Designed to transport Altair/Roomer S. Made of aluminum (plastic basket). Part no. 50346



### 4-Point Powered Spreader Bar

For hospitals and other care environments that use clip slings. Three sizes. Part no. 92441, 92442, 92443



### Lift Charger

Complete with battery and carry case. For charging lifts. Part no. 71470, 71471, 71472, 71473




### Adaptive D-Ring Clip

Allows connection of our portable lifts to existing rail systems. Part no. 53043



### Reacher Arm

Assists when hanging/taking down portable lifts so they do not have to be moved by user. Part no. 50494

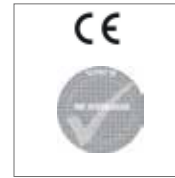

Tested to all relevant standards, including ISO 10535, hoists for the transfer of disabled persons.

## Technical Specifications

Hanger Bars						
Part no.	Description	Length	Weight Capacity	Product Weight	Suggested Use	Suspension
55751	Quick Release Series for Altair/ Roomer S	45cm / 17.7"	300kg / 661lb	1.2kg / 2.6lb	Most lifting situations	2, 3 or 4 <sup>1</sup> -point <sup>1</sup> with 55619
55764	Quick Release Series for Altair/ Roomer S	55cm / 21.7"	300kg / 661lb	1.5kg / 3.2lb	More robust or sensitive to pressure	2, 3 or 4 <sup>1</sup> -point <sup>1</sup> with 55619
72754	Quick Release Series for HeliQ/HeliQ PT	35cm / 13.8"	300kg / 661lb	0.9kg / 2.0lb	Children, Slim Patients	2-point
72755	Quick Release Series for HeliQ/HeliQ PT	45cm / 17.7"	300kg / 661lb	1.1kg / 2.4lb	Most lifting situations	2-point
72757	Quick Release Series for HeliQ/HeliQ PT	60cm / 23.6"	300kg / 661lb	1.3kg / 2.9lb	More robust or sensitive to pressure	2-point
55611	Fixed, for HeliQ/HeliQ PT	35cm / 13.8"	300kg / 661lb	2.0kg / 4.4lb	Children, Slim Patients	2-point
55612	Fixed, for HeliQ/HeliQ PT	55cm / 21.7"	300kg / 661lb	2.0kg / 4.4lb	More robust or sensitive to pressure	2-point
55613	Fixed, for HeliQ/HeliQ PT	60cm / 23.6"	300kg / 661lb	2.0kg / 4.4lb	More robust or sensitive to pressure	2-point
54043	Fixed, for HeliQ/HeliQ PT	45cm / 17.7"	300kg / 661lb	2.0kg / 4.4lb	Most lifting situations	2-point

Hanger Bars for 4-Point Suspension			
Part no.	Dimensions	Weight Capacity	Weight
55619	47cm x 13.4cm / 18.5" x 5.3"	300kg / 661lb	0.8kg / 1.8lb
70800	61cm x 43cm / 24" x 17"	300kg / 661lb	7.2kg / 15.8lb

Accessories			
Part no.	Description	Weight Capacity	Weight
53043	D-Ring	300kg / 661lb	0.1kg / 0.3lb
50494	Reacher Arm	11.2kg/24.6lb	0.5kg / 1.1lb

Continued on next page

This product is manufactured in accordance with all relevant ISO standards.



Follow Us on Social Media



# Technical Specifications

Trolley	
Part no.	50346
Dimensions (L x W)	48.2cm x 40.6cm / 19" x 16"
Height	78.7cm / 31"
Height with Handle	90cm / 35.5"
Wheel Size	12.7cm / 5"
Basket Size (L x W x D)	47cm x 34.3cm x 24cm / 18.5" x 13.5" x 9.5"
Material	Aluminum / Plastic Basket
Product Weight	14.9kg / 32.9lb
Weight Capacity	20kg / 44.1lb
Compatibility	Altair, Roomer S

Lift Charger	
Part no.	Europe: 71470 Great Britain: 71471 Japan: 71472 North America: 71473
Input	100-240V AC, 1.6A Max, 50 - 60Hz
Output	27.8V, 0.8A, DC
Approvals	IP Class IPX0, CE-marked Compliance to EN-60601-1 (Electrical Safety) and EN-60101-2 (EMC)
Protection	Short circuit, polarization, overheating, Class II against electrical shock
Charging Time	6 hours
Components	Cord, Battery Holder, Battery
Connection	Requires connection to a protective earth power socket

Scale (Global)	
Part no.	50590 (Altair/Roomer S) 50591 (HeliQ)
Dimensions (L x W x D)	12cm x 7cm x 16cm / 4.7" x 2.8" x 6.3"
Display	3./cm / 1.5" LCD Panel 5.5 digit display
Battery	6 x 1.5V AAA
Precision Class	2009/23/EC OIML III
Temperature	Ambient: 0°C to +40°C
Approvals	Directive 93/42/EEC, 2009/23/EC CE, FDA D072479/3007074192
Product Weight	1.3kg / 2.9lb
Weight Capacity	300kg / 661.4lb

4-Point Spreader Bar			
	Small	Medium	Large
Part no.	92441 (Manual)	92442 (Electric)	92443 (Electric)
Dimensions (L x W x H)	46.5cm x 45cm x 49.5cm / 18.3" x 17.7" x 19.4"	46.5cm x 59cm x 49.5cm / 18.3" x 23.2" x 19.4"	46.5cm x 73cm x 49.5cm / 18.3" x 28.7" x 19.4"
Input	24V DC, 4A		
Conditions	Ambient Temperature: 10°C to +40°C, Ambient Humidity 30% to 75%, Ambient Pressure 700 hPa to 1060 hPa		
Duty Cycle	-	10%, Maximum 2 minutes ON, 18 minutes OFF	10%, Maximum 2 minutes ON, 18 minutes OFF
Approvals	CE Approved to directive MDD93/42/EEC, In compliance with ISO10535 Models 92442, 92443 IP Rating IPX5		
Product Weight	7.6kg / 16.8lb	10.2kg / 22.5lb	10.6kg / 23.4lb
Weight Capacity	227kg / 500lb		

Scale (North America only)	
Part no.	50348
Dimensions (L x W x D)	16.4cm x 26.5cm x 7cm / 6.5" x 10.4" x 2.7"
Resolution	0 - 200kg = 100g 200 - 320kg = 200g
Battery	4 x LR6/AA
Precision Class	OIML R76 C2/C1.6
Temperature	Ambient: -10°C to +40°C Storage: -40°C to +85°C
Approvals	IP Class: IP54 Directive 90/384/EEC
Product Weight	0.7kg / 1.5lb
Weight Capacity	320kg / 705lb

Upholstered Cover			
Part no.	50866	Material	Polyester
Dimensions (L x W)	36cm x 23cm / 14.2" x 9.1"	Hanger Bar Compatibility	55751, 55764, 72755, 72757
Warnings	Not for use in a wet environment		

This product is manufactured in accordance with all relevant ISO standards.



Learn More About This Product



## Follow Us on Social Media

