

Fresstanding Gantry Lena



Our Lena free-standing gantry system is a perfect alternative when building modification is not possible. It works seamlessly with all our lifts and its flexibility means that it is easily moved to where it's needed most. With a high weight capacity, the Lena provides a safe and easy lifting solution for most patients.

Portable Solution

Sometimes a facility or homeowner does not want to modify the existing building structure to install an overhead rail system or sometimes only a temporary solution is required. That is when the Lena is the perfect option. Portable, easy to use and assemble, it won't affect walls or ceilings. Its compatibility with our Altair series, HeliQ series and Roomer S allow it to be used by a number of patients or lifting situations.

Versatile

The Lena is available in two widths, both of which easily fit over a standard bed. Combined with our overhead lifts and slings, a caregiver can easily lift and move patients without risk of injury. The aesthetic powder

coated white finish ensures the Lena seamlessly fits into any room environment including bathrooms. With optional wheels added, it can easily be transported from room to room.

Quality

Our products are designed and manufactured to meet all the requirements of relevant ISO, UL, IEC and EN standards. At every stage of the development process, every single product is given extensive attention in design, development and manufacturing. This means that a user of this product feels confident with its reliability and the caregiver feels confident with its safety.



Walking Aids



Healthcare Beds



Lifting Solutions



Hygiene



Convertible chairs

Accessories



Rail Rollers

Optional rail rollers with hooks can be used with this product (Part no. 20529 shown). Part no. 20107 with loop included.



Scale

Medical grade scale can be used on all Human Care lifts. Part no. 50590 (Roomer S), 50591 (Altair/HeliQ), 50345 (NA).



Distance Strap

Reinforced nylon thread and metal hooks. Extend lift strap up to 100cm/39.3". Part no. 20124-xx (specify length)



Wheels

Optional wheel kit to assist with portability. Part no. 50872

Lifting Situations



Use to/from a bed

Can be placed over bed for transfers and repositioning with a ceiling lift/sling.




Use to/from a chair

Can be placed over a chair for transfers with a ceiling lift/sling.



Use to/from the floor

Can be placed over a patient in the supine position for transfers/lifting with a ceiling lift/sling.

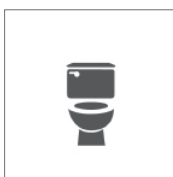

Use to/from the toilet

Can be placed over a toilet for transfers with a ceiling lift/sling.

Features



Powder Coated

Rail: Powder coated extruded aluminum
Poles and feet: Powder coated steel.




Portable Lifts

Compatible for use with our Altair series of ceiling lifts (shown) and our Roomer S.



Stationary Lifts

Compatible for use with our HeliQ series of ceiling lifts.

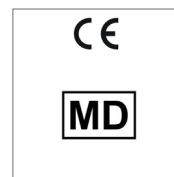

Tested to all relevant standards, including ISO 10535, hoists for the transfer of disabled persons.

Technical Specifications

Part Number	Description	Length	Width	Height	Product Weight	Weight Capacity
71806	Lena Freestand (2.5m)	250cm / 98.4"	145.4cm / 57.2"	229cm / 90.2"	42kg / 92.6lb	300kg / 661lb
71808	Lena Freestand (3.5m)	300cm / 118.1"	145.4cm / 57.2"	229cm / 90.2"	43.5kg / 96lb	300kg / 661lb
50872	Castor kit Lena stand	Includes: 4 pcs (of castors)				

This product is manufactured in accordance with all relevant ISO standards.



Learn More About This Product

