



- Cable-free
- Wide backstrap
- Three colours



Our bestselling classic

neXus 3, is our bestselling rollator with a wide backstrap and robust wheels, making it a perfect everyday-life-companion.



Easy to fold
The best selling neXus series rollators have innovated the way we think about walking aids. With the first-ever patented cross-folding frame, it set the standard for future walking aids design. The neXus 3 rollator features a cross-folding frame, making it exceptionally easy to fold. The soft zippered bag folds together while remaining on the rollator. It can stand independently in a folded position with the bag on the rollator. The cross-folding design also enables the rollator to be partially folded to pass through narrow spaces.

braking system also ensures years of trouble-free use. We have had anatomy and strong visual communication in mind when designing the brake handles. They are comfortable, and the red ball makes the brake handles visible. The cross-folding structure enables superior walking ergonomics with ample walking space. The padded, contoured, and cushioned seat offers a comfortable resting spot. The backstrap can be flipped to allow for two-position seating.



Cable-free and wide backstrap
The cable-free and anatomic brake handles increases safety and comfortability. No cables catch on objects around the home. The cable-free

Available in black, blue and red
The neXus 3 is available in three different sizes; neXus 3 Super Low, neXus 3 Low, and neXus 3 Standard. All sizes are available in Black, Blue and Red (Burgundy) colour.


Features



Padded Backrest
Back comfort while seated. Can be flipped both ways and is included with the rollator.



Foam Cushioned Seat
Offers a comfortable resting spot while walking.

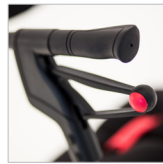


Anatomic Brake Handles
The anatomic brake handles are easy and comfortable to use.

Accessories



Tray
Practical for resting lightweight items. Maximum load 2 kg/4.4lb. Part no. 93316.

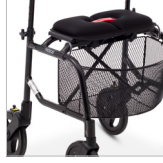

Hard Basket
Basket ideal for shopping. Maximum load 4.5kg/10 lb. Part no. 4510.



Cane or Umbrella Holder
Simplifies transport of a cane or umbrella. Part no. 93317.



Cup Holder
Attaches in seconds without tools. Fits standard size bottle and cups. Part no. 90394.



Curb Climber
Allows the rollator to be slightly lifted onto a curb. Part no. 93315.



Identification Tag Green
Attaches to the rollator in seconds. Part no. 90393.



Soft Zippered Bag
Soft bag that folds with the mobility aid. Maximum capacity 5kg/11lb. Part no. 93280.

Technical specifications

Part Name	Part no.	Handle Height [Min - Max]	Seat Height	Overall Width	Width between Handles	Distance between Rear Wheels	Wheel Size	Product Weight
neXus 3 Super Low	Black 470018	78 cm - 86 cm 30.5" - 34.1"	48 cm 19"	58 cm 23"	45 cm 17.5"	47 cm 18.5"	15 cm x 3 cm 6" x 1.2"	7.6 kg 16.9lb
	Blue 470018B							
	Red [Burgundy] 470018R							
neXus 3 Low	Black 470022	78 cm - 93 cm 30.5" - 36.5"	56 cm 22"	58 cm 23"	45 cm 17.5"	47 cm 18.5"	15 cm x 3 cm 6" x 1.2"	7.7 kg 17.1lb
	Blue 470022B							
	Red [Burgundy] 470022R							
neXus 3 Standard	Black 470024	83 cm - 101 cm 32.5" - 39.8"	61 cm 24"	58 cm 23"	45 cm 17.5"	47 cm 18.5"	15 cm x 3 cm 6" x 1.2"	7.8 kg 17.3lb
	Blue 470024B							
	Red [Burgundy] 470024R							

Maximum User Weight: 200 kg / 441 lb



The neXus 3 is available in three different colours; Black, Blue, and Red [Burgundy].

We improve freedom of movement and quality of care for people.



Human Care USA Inc.
8006 Cameron Street, Suite K
Austin, TX 78754 USA
+1 512 476 7199
info.us@humancaregroup.com

www.humancaregroup.com