



# Alpha Lifting Mat

Foldable  
4 / 6 grips  
Flexible

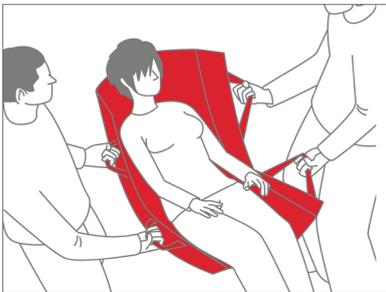

With the Alpha Lifting Mat it is possible to move a client lying in bed towards the top end and on the side, transfer a client, e. g. from a bed to a shower stretcher, and lift up a client who has fallen on the floor.

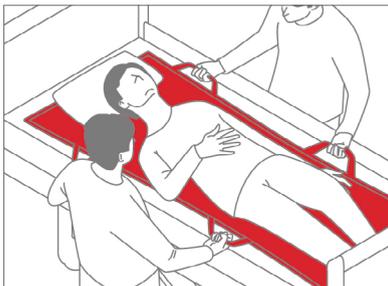

A heavy-duty sliding surface is sewn to the bottom of the transfer sheet. The sliding surface significantly decreases the frictional resistance between the client and the transfer surface. The lateral sections can easily be folded inwards together with the handles, and with the client positioned on his side, the transfer sheet can easily be placed underneath the client. The inner section of the transfer sheet can be removed so that it can be washed or disinfected.

## Technical specification

Product Name	Art. nr.	Dimensions	Material	Max lifting weight
Alpha Lifting Mat with 4 grips	PM-4012	110 x 53 cm 43" x 20,8"	Polyamide, latex-free	200kg / 440lb
Alpha Lifting Mat with 6 grips	PM-4016	150 x 53 cm 59" x 20,8"	Polyamide, latex-free	200kg / 440lb



**PETERMANN**  
part of the Human Care Group

These products are distributed by  
Human Care AB for Petermann GmbH.



Visit [www.humancaregroup.com](http://www.humancaregroup.com) for the latest product information and contact details to national offices and partners.