
neXus 3





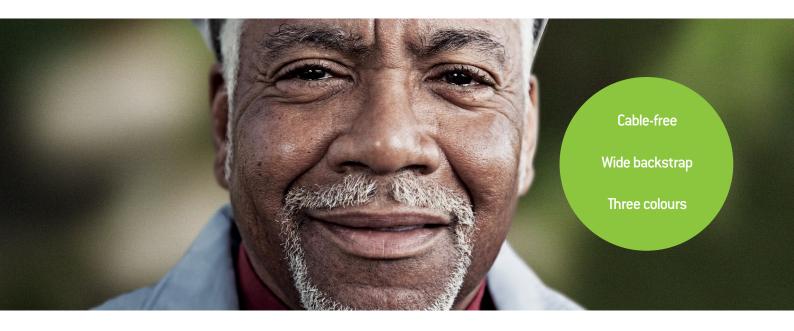

Our bestselling classic

neXus 3, is our bestselling rollator with a wide backstrap and robust wheels, making it a perfect everyday-life-companion.





Easy to fold

The best selling neXus series rollators have innovated the way we think about walking aids. With the first-ever patented crossfolding frame, it set the standard for future walking aids design. The neXus 3 rollator features a cross-folding frame, making it exceptionally easy to fold. The soft zippered bag folds together while remaining on the rollator. It can stand independently in a folded position with the bag on the rollator. The cross-folding design also enables the rollator to be partially folded to pass through narrow spaces. Included with the neXus 3 are the accessories back strap and soft basket.

Cable-free and wide backstrap

The cable-free and anatomic brake handles increases safety and

comfortability. No cables catch on objects around the home. The cable-free braking system also ensures years of trouble-free use. We have had anatomy and strong visual communication in mind when designing the brake handles. They are comfortable, and the red ball makes the brake handles visible. The crossfolding structure enables superior walking ergonomics with ample walking space. The padded, contoured, and cushioned seat offers a comfortable resting spot. The backstrap can be flipped to allow for two-position seating.

Available in black, blue and red

The neXus 3 is available in three different sizes; neXus 3 Super Low, neXus 3 Low, and neXus 3 Standard. All sizes are available in Black, Blue and Red (Burgundy) colour.













Padded Backrest Back comfort while seated. Can be flipped both ways and is included with the rollator.



Accessories

Tray
Practical for resting lightweight items.
Maximum load
2 kg/4.4lb.
Part no. 93316.



Curb ClimperAllows the rollator to be slightly lifted onto a curb. Part no. 93315.



Foam Cushioned Seat Offers a comfortable resting spot while walking.



Hard Basket
Basket ideal for
shopping.
Maximum load
4.5kg/10 lb.
Part no. 4510.



Identification Tag Green Attaches to the rollator in seconds. Part no. 90393.



Anatomic Brake Handles
The anatomic brake
handles are easy and
comfortable to use.



Cane or Umbrella Holder Simplifies transport of a cane or umbrella. Part no. 93317.



Soft Zippered Bag Soft bag that folds with the mobility aid. Maximum capacity 5kg /11lb. Part no. 93280.



Cup HolderAttaches in seconds without tools. Fits standard size bottle and cups. Part no. 90394.

Technical specifications

| Part Name | Part no. | Handle Height (Min - Max) | Seat Height | Overall Width | Width between Handles | Distance between Rear Wheels | Wheel Size | Product Weight |
|---|------------------------------|---------------------------------|--------------|---------------|--------------------------|------------------------------------|---------------------------|------------------|
| | 470018 470018B | 78 cm - 86 cm 30.5" - 34.1" | 48 cm 19" | 58 cm 23" | 45 cm 17.5" | 47 cm 18.5" | 20 cm x 3 cm 8" x 1.2" | 7.6 kg 16.9lb |
| | 470022 470022B 470022R | 78 cm - 93 cm 30.5" - 36.5" | 56 cm 22" | 58 cm 23" | 45 cm 17.5" | 47 cm 18.5" | 20 cm x 3 cm 8" x 1.2" | 7.7 kg 17.1lb |
| neXus 3 Standard Black Blue Red (Burgundy) | 470024 470024B | 83 cm - 101 cm 32.5" - 39.8" | 61 cm 24" | 58 cm 23" | 45 cm 17.5" | 47 cm 18.5" | 20 cm x 3 cm 8" x 1.2" | 7.8 kg 17.3lb |

Maximum User Weight: 200 kg/441 lb







The neXus 3 is available in three different colours; Black, Blue, and Red (Burgundy).



We improve freedom of movement and quality of care for people.