

Fellow Supreme



Suitable for many users, the Fellow Supreme combines high quality ergonomic design with a low centre of gravity for increased safety and stability. The forward facing handles allow the user to walk in a more upright position, which is better for both the back as well as the shoulders.

Unique Design

The Fellow Supreme features a unique fork construction that makes it extremely easy to maneuver, especially in smaller spaces. The flip-up seat allows the user to walk into the seat which allows them closer access to counters while doing work in and around the home.

Ergonomic Grips

The Fellow Supreme's soft braking system combined with the ergonomic design of the handles provide users with reduced hand functions the ability to easily lock the wheels in place and comfortably maneuver the rollator while walking. The Fellow Supreme folds easily with a one-hand operation and

can stand independently while folded. Its vast handle height adjustment range make it suitable for most users.

Practical Accessories

The Fellow Supreme comes with an assortment of accessories to make day-to-day life easier. From a tray to a cane holder, you can easily find an accessory to suit your needs. The grid easily removes for cleaning and maintenance.

HUMAN CARE SWEDEN (HQ)

Årstaängsvägen 21B
117 43 Stockholm
Phone: +46 8 665 35 00
Fax: +46 8 665 35 10
info.se@humancaregroup.com

HUMAN CARE AUSTRALIA

Unit 4, 49 Bertie Street
Port Melbourne, VIC 3207
Phone: +61 1300 132 736
Fax: +61 (03) 9645 4954
info.au@humancaregroup.com

HUMAN CARE CANADA

10-155 Colonnade Road
Ottawa, ON K2E 7K1
Phone: +613.723.6734
Fax: +613.723.1058
info.ca@humancaregroup.com

HUMAN CARE NEDERLAND

Elspeterweg 124
8076 PA Vierhouten
Phone: +31 577 412 171
Fax: +31 577 412 170
info.nl@humancaregroup.com

HUMAN CARE UNITED STATES

8006 Cameron Road, Suite K
Austin, TX 78754
Phone: +512.476.7199
Fax: +512.476.7190
info.us@humancaregroup.com

Accessories/Features



Tray

Practical for resting lightweight items. Maximum load 2kg/4.4lb. Part no. 10667



Soft Bag

Soft carry case ideal for personal items. Part no. 08808



Reflector Kit

Includes 30 self-adhesive reflectors. Part no. 91090



Mesh basket

Ideal for those great shopping adventures. Maximum capacity 4.5kg/10lb. Part no.: 10666



Newspaper Strap

Keeps magazines and newspapers secure while walking. Part no. 10679




Security Lock

Flex cord that protects rollator from theft while unattended. Part no. 91092



Cane holder

Simplifies transport of a cane or umbrella. Part no. 10603

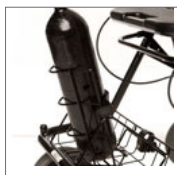

Oxygen Tank Holder

Part no. 10671




One-Hand Brake

Allows control of brakes with one hand. Must be ordered with rollator. Part no. 10669 (right), 10676 (left)



Slow-Down Brake

Prevents the rollator from moving too fast. Part no. 10625 (set), 88121 (right), 88122 (left)

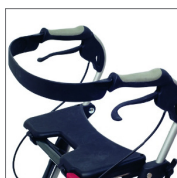

Backrest

Backrest for extra comfort while seated. Part no. 10668

Technical Specifications

Part no.	85400
Handle Height (Min - Max)	73cm - 100cm / 28.7" - 39.4"
Seat Height	59cm / 23.2"
Overall Width	63cm / 24.8"
Overall Length	63.5cm / 25"
Width between Handles	44cm / 17.3"
Distance between Rear Wheels	55cm / 21.7"
Wheel Size	20cm x 4cm / 7.7" x 1.7"
Product Weight	8kg / 17.6lb
Maximum User Weight	150kg / 330.7lb