



The Rebel is a compact rollator ideal for indoor or outdoor use. Its reactive frame is stable, yet easy to maneuver. The Rebel features ergonomic handles that are angle-adjustable allowing the user to find their optimal hand position. With an easy-to-fold system, it can easily be transported.

Perfect Indoors or Outdoors

At home or on the go, the Rebel is the perfect choice. Its reactive frame combined with stable wheels allow the rollator to move smoothly on virtually any surface. Available in an assortment of widths, the Rebel offers a versatile rollator that can be used almost anywhere. The reinforced model for those more robust users features a weight capacity of 250kg/551lb.

Easy Folding System

The Rebel folds easily with a convenient folding handle built right into the seat and stands independently even in this folded position. Its compact size and folding system make it an ideal option for transport and

travel. Shown with the basket and backstrap¹, the Rebel features a full range of accessories to meet the diverse needs of each individual user.

Ergonomic Handles

In order to offer the best ergonomics, the handles feature a stable and comfortable grip that distributes the weight put on them evenly. In addition, the Rebel's handles are angle and height adjustable, so that an optimal hand position can be found for each user.

¹ Backstrap and basket are optional in all countries except for Canada

HUMAN CARE SWEDEN (HQ)

Årstaängsvägen 21B
117 43 Stockholm
Phone:- +46 8 665 35 00
Fax:- +46 8 665 35 10
info.se@humancaregroup.com

HUMAN CARE CANADA

10-155 Colonnade Road
Ottawa, ON K2E 7K1
Phone: +613.723.6734
Fax: +613.723.1058
info.ca@humancaregroup.com

HUMAN CARE NEDERLAND

Elspeterweg 124
8076 PA Vierhouten
Phone: +31 577 412 171
Fax: +31 577 412 170
info.nl@humancaregroup.com

HUMAN CARE UNITED STATES

8006 Cameron Road, Suite K
Austin, TX 78754
Phone: +512.476.7199
Fax: +512.476.7190
info.us@humancaregroup.com

Accessories/Features



Tray

Practical for resting lightweight items. Maximum load 2kg/4.4lb. Part no. 10602



Splash Guard

Protects legs and clothes. Part no. 10608



Reflector Kit

Includes 30 self-adhesive reflectors. Part no. 91090



Slow-Down Brake

Prevents the rollator from moving too fast and enables a controlled walking speed for the user. Part no.: 10620




Light kit

Illuminate your way while walking. Part no. 91091




Security Lock

Flex cord that protects rollator from theft while unattended. Part no. 91092

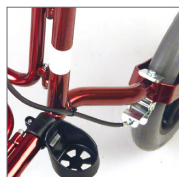

Cane holder

Simplifies transport of a cane or umbrella. Part no. 10603

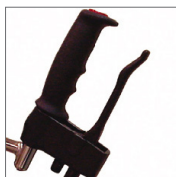


One-Hand Brake

Allows user to control hand functions with one hand (installed left or right). Special order item. Part no. 10652



Bag Kit

Includes one small (5kg capacity) and on large bag (10kg capacity). Part no. 91095 (for 62cm), 91096

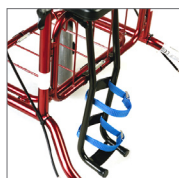

Oxygen Holder

Holds oxygen tube firmly in place. Part no. 10615



Seat Pad

State-of-the-art technology keeps pad cool in summer and warm in winter. Part no. 91093



Backstrap and Basket

Included with rollator in Canada. Backrest Part no. 10621, 10655 (for 72cm) Basket Part no. 10621

Technical Specifications

Description	Part No	Width between Handles	Overall Width	Distance btwn Rear Wheels	Product Weight
Rebel 57	22140 22141 (Low)	37cm / 14.6"	57cm / 22.4"	46cm / 18"	9.5kg / 20.9lb 9.2kg / 20.3lb
Rebel 62	22240 22241 (Low) 71221 (Reinforced)	42cm / 16.5"	62cm / 24.4"	51cm / 20"	9.6kg / 21.2lb 9.4kg / 20.7lb 11.5kg / 25.4lb
Rebel 67	22340 22341 (Low) 71321 (Reinforced)	47cm / 18.5"	67cm / 26.4"	56cm / 22"	9.7kg / 21.4lb 9.5kg / 20.9lb 11.6kg / 25.6lb
Rebel 72	22440 22441 (Low) 71421 (Reinforced)	52cm / 20.5"	72cm / 28.3"	61cm / 24"	9.9kg / 21.8lb 9.7kg / 21.4lb 11.8kg / 26.0lb

Description (All Models)	
Handle Height (Min - Max)	
Regular, Reinforced	75cm - 93cm / 29.5" - 36.6"
Low	68cm - 86cm / 26.8" - 33.8"
Seat Height	
Regular, Reinforced	62cm / 24.4"
Low	54cm / 21.3"
Wheel Size	20cm x 4cm / 7.9" x 1.6"
Seat Size	
Regular, Low	31.8cm x 17.8cm / 12.5" x 7"
Reinforced	40cm x 17.8cm / 15.8" x 7"
Depth folded	25cm / 9.8"
Maximum User Weight	140kg / 330.7lb 250kg / 551lb (Reinforced)

Note: Not all accessories are available for sale in all countries.

This product is manufactured in accordance with all relevant ISO standards.



Follow Us on Social Media



facebook.com/humancareworld



instagram.com/humancaregroup



linkedin.com/company/human-care-ab



pinterest.com/humancaregroup



youtube.com/channel/UC7fyInh3zTjRx0P1IDf2Da0