



The Rebel is a compact rollator ideal for indoor or outdoor use. Its reactive frame is stable, yet easy to maneuver. The Rebel features ergonomic handles that are angle-adjustable allowing the user to find their optimal hand position. With an easy-to-fold system, it can easily be transported.

Perfect Indoors or Outdoors

At home or on the go, the Rebel is the perfect choice. Its reactive frame combined with stable wheels allow the rollator to move smoothly on virtually any surface. Available in an assortment of widths, the Rebel offers a versatile rollator that can be used almost anywhere. The reinforced model for those more robust users features a weight capacity of 250kg/551lb.

Easy Folding System

The Rebel folds easily with a convenient folding handle built right into the seat and stands independently even in this folded position. Its compact size and folding system make it an ideal option for transport and

travel. Shown with the basket and backstrap¹, the Rebel features a full range of accessories to meet the diverse needs of each individual user.

Ergonomic Handles

In order to offer the best ergonomics, the handles feature a stable and comfortable grip that distributes the weight put on them evenly. In addition, the Rebel's handles are angle and height adjustable, so that an optimal hand position can be found for each user.

¹ Backstrap and basket are optional in all countries except for Canada

HUMAN CARE SWEDEN (HQ)

Årstaängsvägen 21B
117 43 Stockholm
Phone:- +46 8 665 35 00
Fax:- +46 8 665 35 10
info.se@humancaregroup.com

HUMAN CARE CANADA

10-155 Colonnade Road
Ottawa, ON K2E 7K1
Phone: +613.723.6734
Fax: +613.723.1058
info.ca@humancaregroup.com

HUMAN CARE NEDERLAND

Elspeterweg 124
8076 PA Vierhouten
Phone: +31 577 412 171
Fax: +31 577 412 170
info.nl@humancaregroup.com

HUMAN CARE UNITED STATES

8006 Cameron Road, Suite K
Austin, TX 78754
Phone: +512.476.7199
Fax: +512.476.7190
info.us@humancaregroup.com

Accessories/Features



Tray

Practical for resting lightweight items. Maximum load 2kg/4.4lb. Part no. 10651



Curb climber

Acts as an aid for walking on pavement and climbing curbs. Part no. 10612



Cross handle

Cross handle with included one-hand brake. Only for Rebel 62. Art.no: left 70844, right 70845



Slow-Down Brake

Prevents the rollator from moving too fast and enables a controlled walking speed for the user. Part no.: 10624




Backstrap and Basket

Included with rollator in Canada. Backrest Part no. 10621, 10655 (for 72cm) Basket Part no. 10601




Security Lock

Flex cord that protects rollator from theft while unattended. Part no. 91092

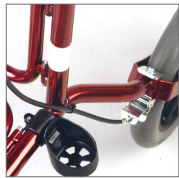

Cane holder

Simplifies transport of a cane or umbrella. Part no. 10603

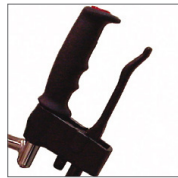


One-Hand Brake

Allows user to control hand functions with one hand (installed left or right). Special order item. Part no. 10652




Bag Kit

Includes one small (5kg capacity) and one large bag (10kg capacity). Part no. 91095 (for 62cm), 91096

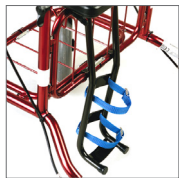

Oxygen Holder

Holds oxygen tube firmly in place. Part no. 10615

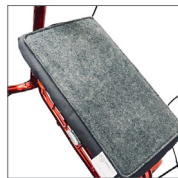

Seat Pad

State-of-the-art technology keeps pad cool in summer and warm in winter. Part no. 91093



Oxygen Holder

Holds larger oxygen tube firmly in place. Part no. 10609

Technical Specifications

	Rebel 67	Rebel 72	Rebel Reinforced 67	Rebel Reinforced 72
Part no.	2234101BR	2244101BR	71321BR	71421BR
Handle Height (Min - Max)	68cm - 86cm/ 26.8" - 33.8"	68cm - 86cm/ 26.8" - 33.8"	75cm - 93cm/ 29.5" - 36.6"	75cm - 93cm/ 29.5" - 36.6"
Seat Height	54cm/21.3"	54cm/21.3"	54cm/21.3"	54cm/21.3"
Overall Width	67cm/26.4"	72cm/28.3"	67cm/26.4"	72cm/28.3"
Overall Length	75cm/29.5"	75cm/29.5"	75cm/29.5"	75cm/29.5"
Width between Handles	47cm/18.5"	52cm/20.5"	47cm/18.5"	52cm/20.5"
Distance between Rear Wheels	52.1cm/20.5"	52.1cm/20.5"	52.1cm/20.5"	52.1cm/20.5"
Wheel Size	20cm x 4cm/ 7.9" x 1.6"	20cm x 4cm/ 7.9" x 1.6"	20cm x 4cm/ 7.9" x 1.6"	20cm x 4cm/ 7.9" x 1.6"
Product Weight	9.5kg/20.9lb	9.7kg/21.4lb	11.6kg/25.6lb	11.8kg/26.0lb
Maximum User Weight	140kg/309lb	140kg/309lb	250kg/551lb	250kg/551lb