

# Orka Shower Trolley

Durable

Easy  
adjustments

Safe



The Orka shower trolley is a robust and comfortable device for showering lying down clients at an ergonomical height for the care giver. The shower trolley is equipped with a hygienic and soft mattress, as well as an excellent pillow and waterproof quality. The Orka is available in junior and bariatric models - offering great aid to a large amount of clients.

## The Orka Shower Trolley Functions

The Orka comes in a hydraulic or electric height adjustment option and includes foldable side rails. The one piece water resistant vinyl covered shower mattress has high sides and a semi-cylindrical head pillow for added comfort and safety. The large castors offer an easy steering system which allows you to lock the front left castor in the direction of steering.

The top end of the stretcher has gas spring assisted incline up to 5°, this prevents pooling around the clients head, as well as facilitates drainage through the integrated mattress drain which includes a flexible hose.

The Orka shower trolley is manufactured from powder coated sturdy galvanised steel to ensure it can withstand the demands of regular wet area exposure. It ensures easy steering, fast adjustments to guarantee flexibility for the caregiver, and a comfortable position for the client.

## Quality

Our products are designed and manufactured to meet all the requirements of relevant standards. At every stage of the development process, every single product is given extensive attention in design, development and manufacturing. This means that a user of this product feels confident with its reliability and the caregiver feels confident with its safety.



Walking Aids



Lifting Solutions



Healthcare Beds



Hygiene



Convertible chairs

# Functions




**Central braking**  
Facilitates maneuvering of the bed in tight spaces.




**Remote control**  
Control, with placing hook, for raising and lowering the shower trolley.



**Drain outlet**  
The drain is positioned in the middle of the mattress's foot end. Plastic or metal option.

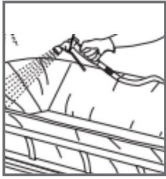

**Comfortable mattress**  
Soft and comfortable mattress with closed corners, that is easy to clean.

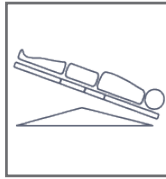

**Trendelenburg**  
16 degrees Trendelenburg / reverse Trendelenburg



**Large castors**  
Large castors makes it easy to manoeuvre



**Head cushion**  
Comfortable head cushion enables the carer to position the clients head safely.



**Head/Side rails**  
Easily unlocked with one hand and not in the way when transferring.



**Drain Hose**  
Comes with a flexible drain hose is 700mm long and can be extended up to 2000mm.

# Technical specifications

Product Name	Width	Length	Safe Working Load	Article Number
Orka Shower Trolley Junior Hydraulic	60 cm	150 cm	170 kg	10068972
Orka Shower Trolley Junior Electrical	60 cm	150 cm	170 kg	10068973
Orka Shower Trolley Standard Hydraulic	65 cm	190 cm	170 kg	10068911
Orka Shower Trolley Standard Electric	65 cm	190 cm	170 kg	10068913
Orka Shower Trolley Bariatric Hydraulic	80 cm	200 cm	250 kg	10068921
Orka Shower Trolley Bariatric Electrical	80 cm	200 cm	250 kg	10068923
Orka Shower Trolley Bariatric Hydraulic	90 cm	200 cm	250 kg	10068931
Orka Shower Trolley Bariatric Electric	90 cm	200 cm	250 kg	10068933



= 170 kg - 250 kg  
= 397 lbs - 550 lbs



*We improve freedom of movement and quality of care for people.*

[www.humancaregroup.com](http://www.humancaregroup.com)



**HUMAN CARE SWEDEN (HQ)**  
Årstaängsvägen 21B  
117 43 Stockholm  
+46 8 665 35 00  
info@humancaregroup.com

**HUMAN CARE CANADA**  
10-155 Colonnade Road  
ON K2E 7K1 Ottawa  
+1 613 723 6734  
info.ca@humancaregroup.com

**HUMAN CARE THE NETHERLANDS**  
Elspeterweg 124  
8076 PA Vierhouten  
+31 577 41 21 71  
info.nl@humancaregroup.com

**HUMAN CARE UNITED STATES**  
8006 Cameron rd, Suite K  
TX 78754 Austin  
+1 5 124 76 7199  
info.us@humancaregroup.com