

Mobi Pro Flexi II



Comfortable
Load capacity 180kg
Ergonomic



The Mobi Pro Flexi II is an ergonomic and compact mobile transfer aid, ideal for people with little leg strength and some trunk balance who have difficulty walking and standing. The Mobi Pro Flexi II is an excellent and safe transfer aid for early mobilisation. It provides support during standing and is an efficient aid for moving patients, improving safety and efficiency.



Comfortable standing aid

The Mobi Pro Flexi II is an efficient and safe standing aid that helps clients to stand up while using their own strength. The height adjustable, padded lower leg/knee supports work to brace the legs of the client during transfer. The sturdy non-slip footboard allows the client to easily position their feet and pull themselves up to a standing position using the crossbar. By simply turning the ergonomic seat supports, the Mobi Pro Flexi II becomes a comfortable chair of which the client can be transferred by. The Mobi Pro Flexi II can effortlessly move a client who weighs up to 180 kilos.

Versatile design

The compact and ergonomic design of the Mobi Pro Flexi II allows a single caregiver to support the transfer from sitting to standing a client to bed, chair or toilet. The Mobi Pro Flexi II is equipped with a mechanical leg spreader with foot pedals, so that the width can easily be adjusted to wide chairs.

Quality

Our products are developed with great care, with the mobility of the care recipient central in a reliable, safe and ergonomic way. All our aids comply with the MDR and all applicable laws, standards and guidelines.



Walking Aids



Healthcare Beds



Lifting Solutions



Hygiene



Convertible chairs

Accessories/Features



Bedside Use

Use to transfer a patient to/from a bed.



Easy Movement

The Mobi Pro Flexi is easy to move by a single caregiver. A small turning circle makes it ideal for smaller spaces.




Padded Knee Support

Comfortable knee pad to provide greater patient comfort and reduce pressure on the knees when standing.



Chair Use

Use to transfer a patient to/from a (wheel)chair.

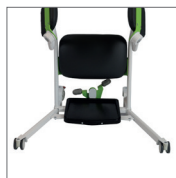


Adjustable Base

Manual leg spread adjustment ensures it fits under most chairs.



Ergonomic Handles

Ergonomic handles for both caregivers and patients.

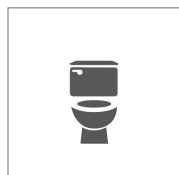

Toilet Use

Use to transfer a patient to/from a toilet.



Padded Seat Supports

Allows the device to be used as a seat, as an alternative to a wheelchair.



Footplate

Durable ABS plastic. Non-slip.

Technical Specifications

Mobi Pro Flexi II			
Article number	100200	Adjustable leg [Min - Max]	49 cm - 92 cm
Width	62 cm	Handle height	101.5 cm
Length	90.5 cm	Weight	35.5 kg
Height	102 cm	Max. Lift capacity	180 kg
Diameter of wheels with brake	10 cm	Diameter of wheels without brake	7.5 cm



We improve freedom of movement and quality of care for people.

www.humancaregroup.com



HUMAN CARE SWEDEN (HQ)

Årstaängsvägen 21B
117 43 Stockholm
+46 8 665 35 00
info@humancaregroup.com

HUMAN CARE CANADA

10-155 Colonnade Road
ON K2E 7K1 Ottawa
+1 613 723 6734
info.ca@humancaregroup.com

HUMAN CARE THE NETHERLANDS

Elspeterweg 124
8076 PA Vierhouten
+31 577 41 21 71
info.nl@humancaregroup.com

HUMAN CARE UNITED STATES

8006 Cameron rd, Suite K
TX 78754 Austin
+1 512 476 7199
info.us@humancaregroup.com